




Caring Vision[®]

Innovative Mammography Solutions for 20 years

The new Rolloscope MD[®] model accommodates digital LCD monitors and prior film comparisons in an ergonomic viewing system.

With the Rolloscope's "One Touch System" the backlight to the monitors is controlled and minimized. Push just one button and the multiviewer's lights automatically dim, the belt advances, the shutters mask to the new frame's film configuration and the unit brightens in about one second.

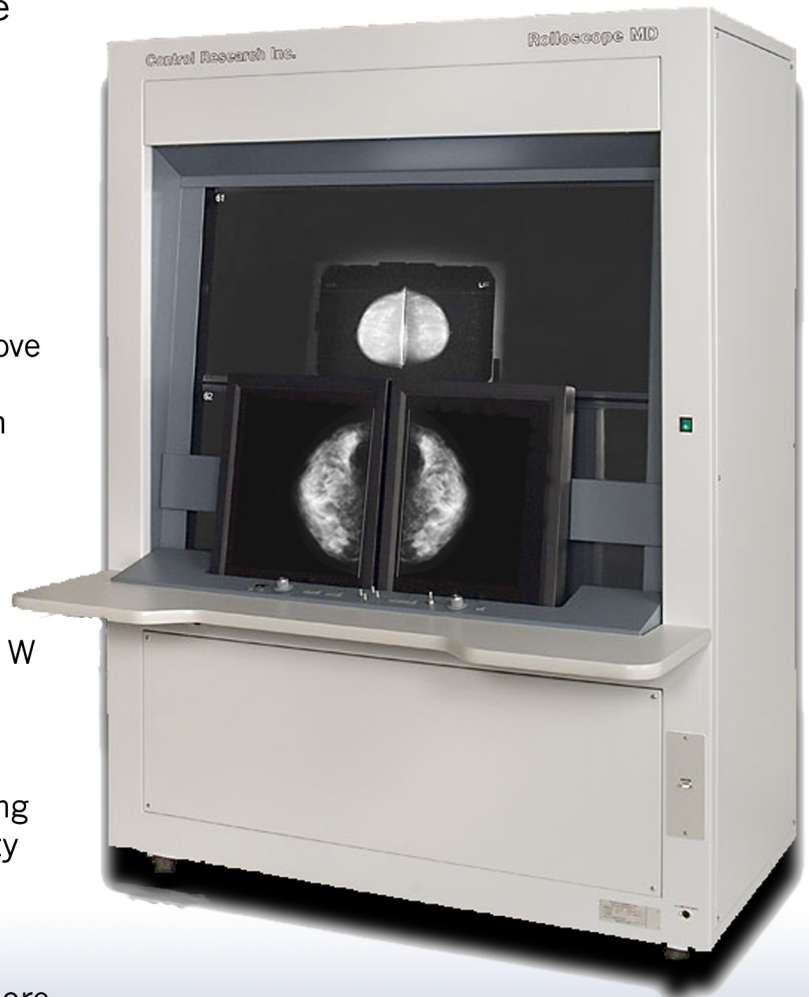

The Rolloscope MD has 101 frames above the LCD monitors. These frames accommodate 6-over-6 of the 18x24cm films and 4 across of the 24x30cm mammos or 3 across of 14x17" MRI or large ultrasound films above monitors.



Without monitors it would be twice the capacity above. The unit measures 56" W x 78.75" H x 30" D (the shelf adds an additional 15" depth at the corners).



With its superior luminance and masking technology, ease of use and high quality construction, the Rolloscope MD offers exceptional value and performance. Please visit our website, www.broadwest.com, for a video and more information.

T. (303)722-1500
(800)232-2948
F. (303)722-1507
info@broadwest.com
www.broadwest.com