

Classic Magnifiers

Round 5 inch diameter magnifiers, all with 2x magnification. Optical acrylic or crown glass lenses.

RM-5A ACRYLIC lens with 4 power inset & soft ergonomic grip, 9 oz.

New **RM-5G** GLASS lens with soft ergonomic handle, extra wide view, 15 oz.

RM-5GW GLASS lens with solid walnut handle & brass frame, 16 oz.



MammoMag™ Viewers

MammoMag enhances perception of subtle radiographic detail. By blocking ambient light, contrast perception is improved and fine glandular tissue differences are observed more readily. It magnifies and creates a focused field of view. Choose from three interchangeable glass binocular lenses. Dimensions: 9.5" x 6" x 3", 11 oz.

MS2 1.5x magnification, 20" focal range.

MS3 1.75x magnification, 14" focal range.

MS4 2x magnification, 10" focal range (body of MammoMag is 9.5").

Glass replacement lenses: **GRL2** (1.5x), **3** (1.75x) or **4** (2x).



Aspheric Dome Magnifier

Solid acrylic dome prefocuses and magnifies the image by a factor of 4. The base of the dome is slightly concave to resist scratches. Ideal for calcifications, lesions, subtle tissue distortions, film quality control and slides. A felt lined wooden base and flannel bag are provided.

4X-D 2.5" diameter, 5oz.

4X-DT 2.5" tethered, 5.5oz., affixes to viewbox.

4X-D3.5 3.5" diameter, 13oz.

New 30mm detachable clear scale included with each dome.



Headlamp Magnifier

The adjustable visor headlamp has an optical glass magnifying lens and a 8,000 candlepower xenon light. Useful for ductography, sialograms or other interventional procedures. Padded leather band for comfort and friction locks for secure positioning. Operates on two AA batteries which are included. Three magnification levels available:

LMV-2 1.5x, 20" focal length **LMV-3** 1.75x, 14" focal length **LMV-4** 2x, 10" focal length, 8oz.

Also available without light: **MV-2** 1.5x, **MV-3** 1.75x, **MV-4** 2x, 5oz.



Bright-Eye™ Spot with Iris

Bright-Eye spot with adjustable iris diaphragm allows a 1/4" to 4" diameter field of light. Excludes extraneous light to easily pinpoint area of concern on the film. Foot-operated dimmer provides variable light intensity. Rotates on base and tilts back and forth. 110 watt reflector lamp with a rough service filament provided for high luminance and long service.

Model **SP-2**.



Technologist's Viewer


Sturdy & economical illuminators for checking films in mammo/biopsy room or near processor. Two light intensities, 4,000 nit luminance, 14watt lamps, positive film grips. Holds 18x24cm or 24x30cm mammo format.

T-2 Single bank: 2 films of 24x30cm, 5 lamps
Viewing surface: 48 x 30cm.

Dimensions: 16 3/8"H x 19 1/8"L x 4"D.

T-4 Double bank: 4 films of 24x30cm, 10 lamps. Viewing surface: 96x30cm.

Dimensions: 16 3/8"H x 39 1/4"L x 4"D.

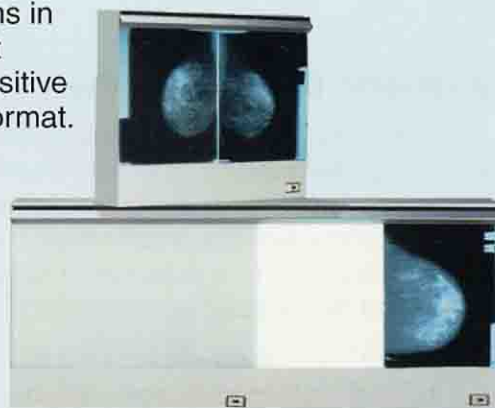

Distributed by:

BROADWEST CORPORATION

304 Elati Street Denver CO 80223

Tel 303-722-1500 Fax 303-722-1507

8 0 0 - 2 3 2 - 2 9 4 8

www.broadwest.com