

# Carestream

## Image Suite Software V2

with DIRECTVIEW Vita, DIRECTVIEW Classic  
and Point-of-Care 140 and 360 CR Systems

## SIMPLE. AFFORDABLE. POWERFUL.

The perfect combination for workflow productivity – a complete digital image acquisition and management solution designed specifically for smaller facilities. This "mini-PACS" software solution will enable digital imaging in your facility and streamline your workflow. Combined with a robust and compact CR system, you'll get the high quality diagnostic images you need that fits perfectly in your office space – all at an affordable price.



## SIMPLE. COMPACT.

# IMAGE SUITE V2 IS THE PERFECT SOLUTION.

## Simple - Plug it in and go

Quickly and easily acquire, process, edit and send images for storage or printing. It's simple to install and training can be completed in less than a day.

## Affordable

Created especially for private practices and smaller clinics to offer an economical PACS solution.


## Powerful

All in one system with superior automatic image processing.

## Flexible

Easily customized to meet your specific business needs. Get the features you need now, and add capabilities in the future as your practice grows.

# PERFECT COMPANIONS – CARESTREAM'S IMAGE

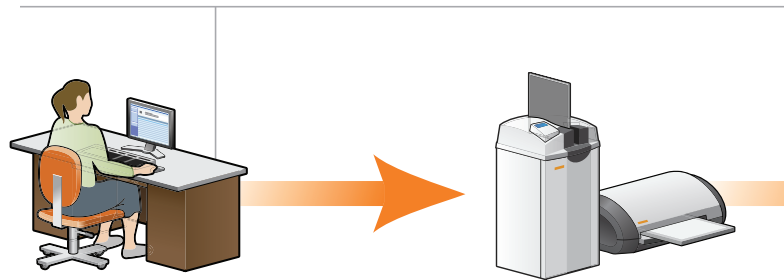

## SCHEDULE

### SCHEDULE

#### Speed up the Registration Process - From Anywhere

- Register patients remotely with Web-based software
- Request multiple images in a single order
- Flexibility with body parts & views – add new, edit existing or apply custom processing
- Connect with DICOM or HL7 Compliant EMR/RIS systems for simplified workflow

### ACQUIRE

#### Perfect Pair with DIRECTVIEW CR Systems for Fast, High Quality Images

- Supports CARESTREAM portfolio of CR Scanners\*
- Includes Enhanced Visualization-Automatic (EV-A) image processing for consistent high quality images every time
- Supports Mammography with DIRECTVIEW Classic CR (outside US)
- Supports Long-Length Imaging\*\*

DIRECTVIEW  
Vita CR



# SUITE SOFTWARE + CR SCANNERS + LASER PRINTERS



## VIEW/EDIT

### Easy to Adjust for Top Quality Images

- Supports eight concurrent web viewer sessions for increased flexibility
- Image processing can be tuned with the click of a button, using simple slide bars to adjust brightness, contrast or detail
- Delivers full bit-depth images at the host workstation
- Advanced measurement tools for chiropractic review
- Tools for mammography review (outside US)

## DIAGNOSE/REPORT

### Faster Diagnoses, Better Patient Care

- Dual monitor support on Web clients enables increased efficiency by allowing the viewing of multiple images at once
- Process full bit-depth (raw/full-fidelity) images at the host workstation, including ability to control brightness, latitude, contrast, sharpness and noise
- Receive images from CT, MR, ultrasound and additional CRs with the optional DICOM store SCP software

## PRINT/ARCHIVE

### Share and Store Images Efficiently

- Flexible security with on-demand or scheduled backup
- Storage of both “for processing & for presentation” images gives you the most flexibility for reprocessing stored images – both simple and advanced adjustments – unlike any other PACS
- Perform backups to CDs and removable USBs on demand for off-line storage
- Usability of images for remote reporting with JPEG or DICOM export tools
- Support multi-format print composition with independent image editing capability

\* DIRECTVIEW Classic CR, Vita/Vita SE/Vita LE or PoC 140/360  
\*\* Long-Length Imaging available with Classic or Vita/Vita SE/ Vita LE


**DRYVIEW**  
**5700 Laser Imager**  
Exceptional affordability, great simplicity.



# THE IMAGE SUITE SOLUTION

A complete digital image acquisition and management solution.

DIRECTVIEW  
Classic CR

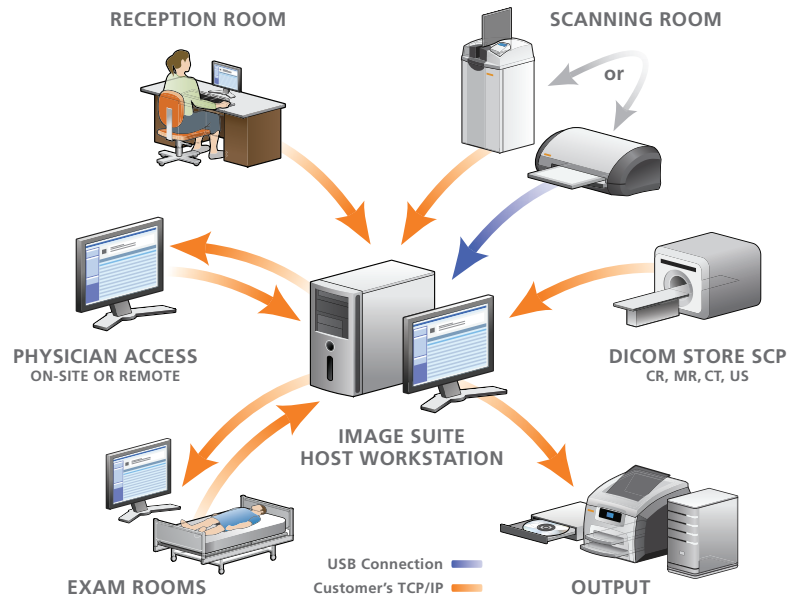

## A Partner That Cares as Much as You Do.

Carestream's comprehensive service offerings, CARESTREAM Smart Services, give you peace of mind over the lifecycle of the equipment. From planning, installation, preventative maintenance, training and service, our certified Service Solution Specialists keep your equipment running at peak performance so you can focus on what is most important – your patients.

## Available Options for Image Suite Software V2

- Clinical Report Software
- Web Viewer License
- Long-Length Imaging Software
- DICOM Store Service Class Provider (SCP)
- Advanced Measurement Tools Set
- Mammography Feature Software (Outside US only)
- Veterinary Application

## Minimum PC Requirements

- Windows 7 Professional Operating System, 32 bit
- Intel Pentium 4 Processor @ 2.8GHz
- 4GB DRAM memory
- 300GB hard disk – must be partitioned for 2 storage directories
- Ethernet card 100 MBp/s (2 each)
- CD/DVD writer
- 17" display of 1024 x 1280 resolution

Distributed by:  
CAN-med Healthcare

East 1 866 277 3577  
West 1 800 661 5671

[www.canmedhealthcare.com](http://www.canmedhealthcare.com)

