


Move forward with the leader in point-of-care carts

INTRODUCING THE NEXT GENERATION!

StyleView[®] Carts

Designed for Nurses.
Engineered for IT.

 More motion. Less effort. **Feel the difference.**



You asked. We listened. We delivered—a point-of-care cart that provides even more caregiver efficiency, patient safety and quality care. Ergotron, the proven leader in point-of-care computing, has refined the StyleView Cart. Built with the same unparalleled features of the original StyleView, the latest version lifts point-of-care computing to a whole new level.

1. **NEW!** Five-inch (13 cm), independent monitor height adjustment increases comfort and productivity for personalized viewing, including for those wearing progressive lenses. Rotates from portrait to landscape mode. Compatible with Tablet PCs (LCD carts only)
2. 20-inch (51 cm) height adjustment for comfortable data entry, sitting or standing
3. Constant Force (CF™) patented motion technology provides fluid, one-touch height adjustment
4. **NEW!** Secure, locking medication drawer controlled by electronic keypad—enhancing caregiver efficiency and patient safety. Removable cassette is configured with six bins (each with a clear lid and ID label)
5. Barcode scanner holder provides easy access to hand-held scanner
6. **NEW!** Wrists remain in a neutral position during data entry with negative-tilt keyboard tray—avoid repetitive-stress injuries and increase comfort and productivity
7. **NEW!** Expansive cable management storage area
8. **NEW!** Easier mousing and typing with rigid keyboard platform
9. **NEW!** Long-life power performance supported with an industry-leading two-year warranty* on power system. Five-year warranty on cart
10. **NEW!** Upgraded StyleLink™ software improves power system management tools, helping extend the life of the battery (reducing costs) by offering battery cycle count, email alerts and free software updates. Windows Vista™ compatible (powered carts only)
11. **NEW!** Easy-to-clean external 10-foot (3 m) coiled cord included with powered carts in North America and China. Powered carts elsewhere include a convenient retractable cord reel
12. Easy-to-use cart comes fully assembled

* Includes all electrical components except battery

Standard Cart Features Worksurface (except Dual Display Cart), basket and rear handle, scanner holder, cable management system, universal CPU holder (LCD models only), full-size keyboard tray with left/right mouse tray, gel wrist rest, mouse holder, height position lock, high-performance casters with dual wheels for carpet or hard floor use. Powered carts include 55Ah, 12VDC, Sealed Lead Acid (SLA), Absorbed Glass Mat (AGM), Deep Cycle, Valve Regulated type battery.



LCD Cart powered with med drawer



Notebook Cart powered with med drawer

Notebook Cart non-powered with med drawer

Notebook Cart non-powered


LCD Cart non-powered

Dual Display Cart powered

Dual Display Cart non-powered



Smart Power using StyleLink Software for StyleView Carts

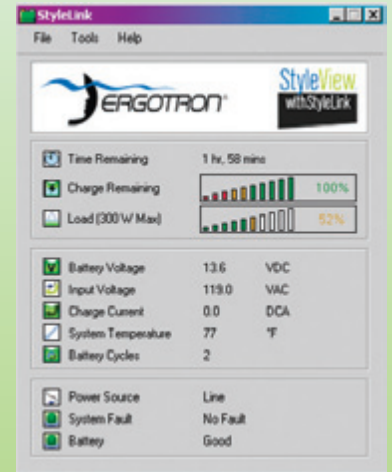

Reduce downtime
Maximize productivity
Reduce cost
Maximize ROI

EXTEND ASSET LIFE CYCLE

Free, downloadable software upgrades for ongoing battery-life optimization

Extends battery life by managing usage based on runtime performance

PROACTIVE POWER SYSTEM MANAGEMENT

Eliminate surprises through use of diagnostic tools, simplifying power system management and reducing time required to fix problems

Email alerts sent by each cart at different stages of battery usage facilitate proactive power system management

StyleLink is included with all powered carts and is Windows Vista compatible. Free online upgrades ensure you'll be ready for future hardware upgrades

REDUCE TOTAL COST OF OWNERSHIP

Reduce costs for maintenance; eliminate need for expensive onsite technical support by leveraging diagnostic capabilities. Simply email StyleLink's report to Ergotron Tech Support and they will help troubleshoot power system performance issues



The power system has a standard two-year warranty which leapfrogs over warranties offered with competitive carts. Ergotron's upgraded StyleLink software enhances power system management, helping extend battery life. For patient and caregiver safety, the cart and power system are UL and EN 60601-1 compliant. A programmable, auto-locking keypad medication drawer is included, making StyleView both a point-of-care cart and a medication-delivery system. StyleView bridges point-of-care computing and medication delivery, reducing healthcare costs and improving patient care.

Improve viewing comfort and enhance productivity through StyleView's expansive range of effortless height adjustment. Twenty inches (51 cm) of sit-to-stand motion is combined with an additional five inches (13 cm) of independent monitor height adjustment and 25° screen tilt, which benefits all users, including those who wear progressive lenses. The keyboard tray provides up to 12° of negative tilt—recommended by ergonomists for the proper typing position. This reduces the risk of repetitive stress injuries. An enclosed storage area makes it easy to manage cables.

	Cart	Part Number	VAC	Hz
Global				
Non-powered	LCD	SV21-91075	—	—
	Notebook	SV21-91007	—	—
	Notebook w/med drawer	SV21-91008	—	—
	Dual Display	SV21-31645	—	—
North America				
Canada	LCD	SV22-92023	120	60
	Notebook	SV22-92026	120	60
	Dual Display	SV22-22009	120	60
USA	LCD	SV22-92023	120	60
	Notebook	SV22-92026	120	60
	Dual Display	SV22-22009	120	60
European Union				
Europe*	LCD	SV22-92014	230	50
	Notebook	SV22-92011	230	50
	Dual Display	SV22-22017	230	50
Ireland United Kingdom	LCD	SV22-92015	230	50
	Notebook	SV22-92012	230	50
Asia Pacific	LCD	SV22-92015	230	50
	Notebook	SV22-92012	230	50
	Dual Display	SV22-22018	230	50
	Dual Display	SV22-22018	230	50
Hong Kong Malaysia Singapore	LCD	SV22-92024	230	50
	Notebook	SV22-92025	230	50
China	LCD	SV22-92016	230	50
	Notebook	SV22-92013	230	50
Australia New Zealand	LCD	SV22-92016	230	50
	Notebook	SV22-22019	230	50



LCD Size	LCD Weight	VESA
≤ 22.9" ≤ 58 cm	6–18 lbs 2.7–8.2 kg	MIS-D

Notebook Size*	Notebook Weight	VESA
≤ 22.9" ≤ 58 cm	5–12 lbs 2.3–5.4 kg	N/A

Dual Display Size	Dual Display Weight†	VESA
≤ 15.4" ≤ 39 cm	≤ 17 lbs ≤ 7.7 kg	MIS-D

Tablet PC Size	Tablet PC Weight	VESA
≤ 22.9" ≤ 58 cm	5–12 lbs 2.3–5.4 kg	MIS-D

*Maximum physical dimensions of a notebook used in the cart 17" x 12.5" x 1.5" (43.1 x 31.7 x 3.8 cm)
† Per LCD pivot

* Contact Ergotron Sales Representative for country-specific configuration

All the details covered

StyleView Accessories

These optional accessories, including lockable medication drawers and utility shelf, are engineered to increase productivity and reduce errors.



Carpet Casters
97-254



Notebook Drawer
68-122-201



Sharps Container and Bracket
60-570-201



Storage Basket
99-160-063



Gel Wrist Rest
85-043



IV Pole
60-566-003

Interface Brackets for Motion Computing Tablet PCs



**Interface Bracket
for Motion Computing
FlexDock (FlexDock not
included)**
60-499-200 black

Compatible with any
LCD mount that uses a
VESA MIS-D mounting
hole pattern



**Interface Bracket
for Motion Computing
LE1600 Tablet PC**
60-579-200 black

Compatible with any
LCD mount that uses a
VESA MIS-D mounting
hole pattern




Locking Drawers
68-124-201



Utility Shelf
50-192-201



T-Slot Bracket
60-575-003

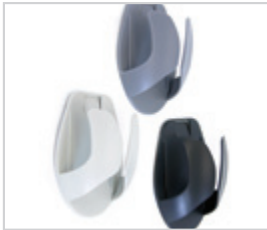

Mouse Holder
99-033-099 grey
99-033-064 dark grey
99-033-085 black



External Cord Reel 10' cord,
110V 60 Hz
83-065



**LX Notebook Tray
with Pivot Extension**
28-510-200



Tablet PC Pivot Extension
47-095

COMING SOON!

5" (13 cm) Monitor Lift Upgrade Kit
99-167-057

Interface Bracket for Motion Computing C5 Tablet PC
60-588-200

Ergotron understands the human health and environmental risks associated with substances covered by the Restriction of Hazardous Substances (RoHS) directive and has eliminated the use of these substances in its global product lines

© 2007 Ergotron, Inc. rev. 04/19/2007 870-05-072-EA Content is subject to change without notification
Some products may not be available in some regions of the globe—visit our website for more information

Americas Sales and Corporate Headquarters

St. Paul, MN USA
(800) 888-8458
+1-651-681-7600
www.ergotron.com
sales@ergotron.com

EMEA Sales

Amersfoort, The Netherlands
+31 33 45 45 600
www.ergotron.com
info.eu@ergotron.com

APAC Sales

Singapore
www.ergotron.com
info.apac@ergotron.com

Worldwide OEM Sales

www.ergotron.com
info.oem@ergotron.com

