


# LIFESAVER<sup>®</sup> Resuscitation Bags

*Manual Labor Never Felt This Good*

The LIFESAVER single patient use manual resuscitator's non-slip grip and quick rebound provide the performance and reliability you need to get a handle on emergency resuscitation.

With one-piece construction to eliminate misassembly and disconnects, a tamper-evident storage bag and numerous configurations to choose from, you'll feel good about the LIFESAVER resuscitation bag from Hudson RCI.

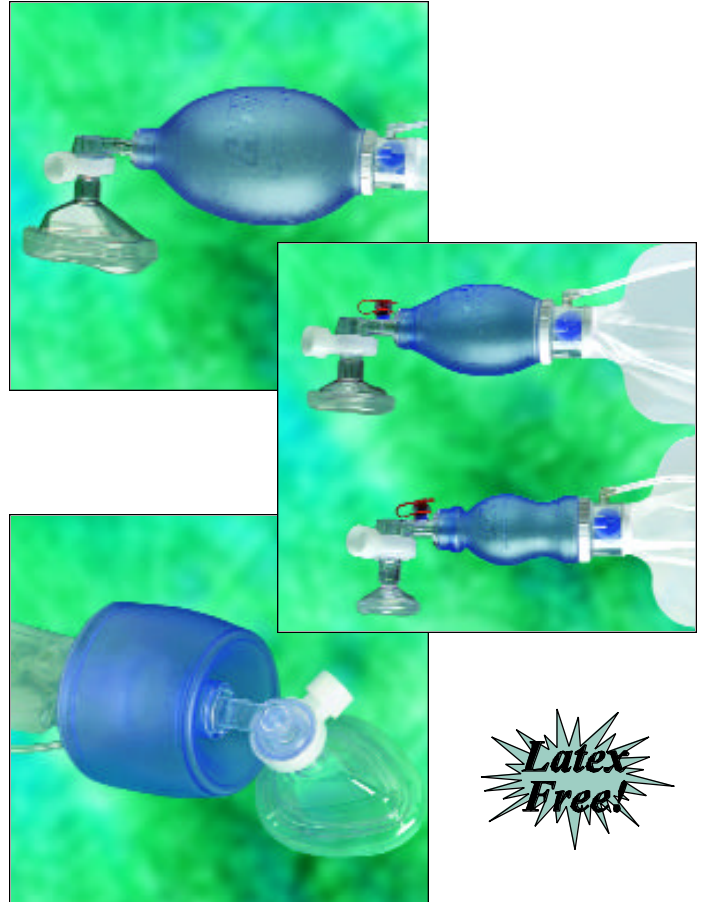

## Features

- Swivel patient connector minimizes trauma to airway and allows easy position changes
- Unique right-angle oxygen tubing connector swivels 360 degrees to facilitate positioning and minimize tubing tangle
- Integral 30mm OD flow diverter allows connection to optional, adjustable PEEP valve
- Foldable adult bag collapses for convenient, compact storage
- Designed to meet all ISO and ASTM resuscitation bag performance standards
- Pediatric and infant bags include capped pressure monitoring port and 40cm H<sub>2</sub>O pressure relief pop-off
- Adult bag available with optional 60cm H<sub>2</sub>O pressure relief and valved mask
- Individually packaged in tamper-evident drawstring bag

## Specifications

Patient Connector	15mm I.D./22 mm. O.D.	<b>Optional PEEP Valve</b>	
Expiratory Connector	30mm. O.D.	• Connector	30mm I.D.
Operating Temperature	-18°C to 50°C.	Range	5 to 20 cmH <sub>2</sub> O ± 2.5 cmH <sub>2</sub> O
Storage Temperature	-40°C to 60°C	• Fits all LIFESAVER bags	
	<b>Infant</b>	<b>Pediatric</b>	<b>Adult</b>
Bag Volume (approx.)	300mL	500mL	1800mL
Stroke Volume	up to 35 mL	up to 300 mL	up to 600mL
Tidal Volume	20mL@60 breaths/min.	300mL@20 breaths/min.	600mL@20 breaths/min.
Dead Space (approx.)			
• Connector	7mL	7mL	7mL
• Mask	18mL	40mL	120mL
Inspiratory Resistance			
• 5 LPM	<1.1cmH <sub>2</sub> O/L/sec	<1.1cmH <sub>2</sub> O/L/sec	n/a
• 50 LPM	n/a	<3.5cmH <sub>2</sub> O/L/sec	<3.5cmH <sub>2</sub> O/L/sec
Expiratory Resistance			
• 5 LPM	<3.0cmH <sub>2</sub> O/L/sec	<3.0cmH <sub>2</sub> O/L/sec	n/a
• 50 LPM	n/a	<3.7cmH <sub>2</sub> O/L/sec	<3.7cmH <sub>2</sub> O/L/sec
Pressure Delivery Limit	40cm H <sub>2</sub> O±5cm H <sub>2</sub> O	40cmH <sub>2</sub> O±5cmH <sub>2</sub> O	60cmH <sub>2</sub> O±10cmH <sub>2</sub> O

Testing performed in accordance with ISO8382: 1988(E), Section A.5.7

## Ordering Information

Catalog No.	Size/Description	Mask Included	Peep Valve	Pressure Pop-Off	Case Qty.
5361	Infant Bag	n/a		40cmH <sub>2</sub> O	6
5362	Infant Bag	Cat. No. 1271		40cmH <sub>2</sub> O	6
5364	Infant Bag	Cat. No. 1271	•	40cmH <sub>2</sub> O	6
5366	Pediatric Bag	n/a		40cmH <sub>2</sub> O	6
5367	Pediatric Bag	Cat. No. 1273		40cmH <sub>2</sub> O	6
5369	Pediatric Bag	Cat. No. 1273	•	40cmH <sub>2</sub> O	6
5371	Adult Bag	n/a		n/a	6
5372	Adult Bag	Cat. No. 1276		n/a	6
5374	Adult Bag	Cat. No. 1282	•	n/a	6
5377	Adult Bag	Cat. No. 1276		60cmH <sub>2</sub> O	6
5380	Adult Bag	Cat. No. 1282		n/a	6
5387	Adult Folded Bag	Cat. No. 1276		n/a	6
5385	PEEP Valve 30mm				10
5383	PEEP Valve 22mm				10



27711 Diaz Road • P.O. Box 9020 • Temecula, CA 92589-9020 USA  
 Tel. 800-848-3766 • 951-676-5611 • Fax 951-676-1578  
 www.hudsonrci.com