

QUICKVUE In-Line® Strep A TEST

Test and Treat Today!™

Strep throat affects all age groups, but is most common in children between the ages of 5 and 15 years of age. Physicians diagnose strep throat based on symptoms, physical findings and culture.

QuickVue In-Line® Strep A is a truly one-step test, which allows for the rapid detection of group A streptococcal antigen, directly from patient throat swab specimens.

QuickVue In-Line Strep A's unique In-Line extraction method reduces hands-on time, and provides quick results allowing testing to take place right at the point of care. Extraction begins instantaneously and the two color test result gives clear readability.



Rapid Diagnosis and Treatment

- Results in 5 minutes, or less
- Most patients tested can be diagnosed and treated in just one visit.

Excellent Performance*

- 92% sensitivity
- 99% specificity
- 96% overall agreement with SBA culture

Easy to Use

- One easy step via In-Line extraction.
- Two-color endpoint for easy to read and interpret results.
- Single reagent method allows for less chance of operator error.

Convenient Packaging

- All components necessary to perform a test are contained in the kit, including swabs.
- Room temperature storage (14 months from date of manufacture).
- Built-in Positive and Negative controls confirm the test has performed properly.
- Positive and Negative external controls provided in the kit.
- Individually packaged for convenient single-use testing and disposal.

*In-house study conducted at Quidel. Refer to Package Insert for additional performance claims.


QUIDEL®
CORPORATION

For more information on QuickVue In-Line Strep A test kits, visit our website or call your Quidel Account Manager at 1.800.874.1517 (U.S. only), or 858.552.1100 (outside U.S.).


Quidel Corporation • Worldwide Headquarters • 10165 McKellar Court • San Diego, CA 92121 USA • www.quidel.com

Procedure:

1
MUST use the swabs provided in the kit

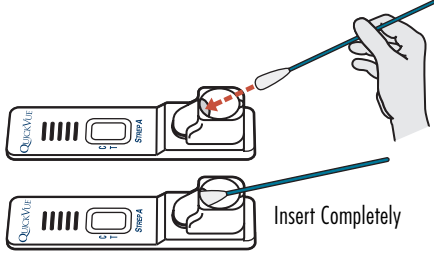


2

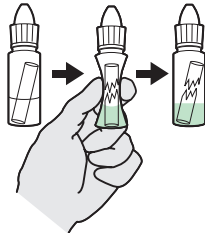


Squeeze **ONCE** to break.

3



Shake

4

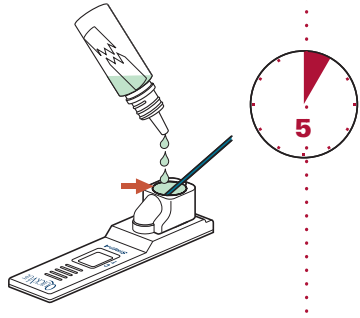

Quickly Fill to Rim



To order: Call your Quidel Account Manager at 1.800.874.1517

QuickVue In-Line® Strep A – 25 Test Kit:
Catalog #00343



For more information on QuickVue In-Line Strep A test kits, visit our website or call your Quidel Account Manager at 1.800.874.1517 (U.S. only), or 858.552.1100 (outside U.S.).

Quidel Corporation • Worldwide Headquarters • 10165 McKellar Court • San Diego, CA 92121 USA • www.quidel.com

Distributor Information: