



## > LIGHTING

### Exam Lights

> Exam Light III	136
> LS-135	136
> LS-150	137
> Professional PenLite	137

### Skin Surface Microscope

> 3.5 V EpiScope Skin Surface Microscope	138
--	-----

### Minor Procedure Lighting

> LS-200	139
----------	-----

### Diagnostic Headlights

> Solid State Portable Headlights	140
-----------------------------------	-----

### Minor Procedure Headlights

> Solid State Procedure Headlights	141
> Solarc MFI Headlight	142

## LIGHTING



Durable Fiber Optic Light Pipe

Sealed Focusing Sleeve



**48859**  
Universal Table Mount



**48760**  
Neonate Fiber Optic Examination Light III



**48605**  
Focusing Sleeve

**48700**  
Exam Light III Fiber Optic Light, Complete shown with **48850** Mobile Stand

## Exam Light III™

Flexible fiber optic lighting that lets you take brighter, whiter light anywhere without the heat.

- > Mounts easily on a wall, stand, or examination table
- > Durable fiber optic light pipe bends to almost any position and transmits cool illumination
- > Easily positioned small head leaves clear view of task areas
- > Sealed focusing sleeve keeps optics free of dust
- > Provides homogeneous spot and variable size
- > Ideal for gynecological examinations and minor procedures
- > Disposable sheaths minimize cross-contamination

- 48700** Exam Light III Fiber Optic Light, Complete
- 48780** Fiber Optic Cable Exam Light
- 48760** Neonate Fiber Optic Examination Light III (complete with infrared filter; no focusing sleeve required)

## Exam Light III™ Accessories

<b>48850</b> Mobile Stand for Exam Light III	<b>04200-U</b> Halogen Lamp
<b>48859</b> Universal Table Mount	<b>52640</b> Disposable Sheaths for Exam Light (25/box, 5 boxes/case)
<b>48200</b> Fiber Optic Light Pipe	<b>52641</b> Disposable Sheath Dispenser
<b>48210</b> Fiber Optic Neonate Light Pipe	<b>48220</b> 7.5 ft Flexible Fiber Optic Cable with Focusing Sleeve
<b>48605</b> Focusing Sleeve	



Molded Lexan Head

Long-Lasting, 35-Watt Halogen Lamp

Dual-Intensity Control Switch



**44101**  
Table Clamp Accessory for 44310



**44310**  
LS-135 Examination Light, Wall/Table Mount Model



**44350**  
Caster Base Accessory for 44300

**44300**  
LS-135 Examination Light, Floor Stand Model

## LS-135™

A unique blend of intense light output, mechanical stability, and affordability.

- > Long-lasting, 35-watt halogen lamp produces significantly whiter, brighter light than ordinary goose neck units
- > Dual-intensity control switch allows instant adjustment of light output
- > Molded Lexan® head does not conduct heat, keeps luminaire cool to the touch
- > Encased, flexible 24 inch arm enables one-hand positioning—not possible with standard goose neck units
- > Multiple installation options to fit your needs

- 44300** LS-135 Examination Light, Floor Stand Model
- 44310** LS-135 Examination Light, Wall/Table Mount Model

## LS-135 Accessories

- 44350** Caster Base Accessory for 44300
- 44101** Table Clamp Accessory for 44310
- 04430-U** 35 W Halogen Replacement Lamp for LS-135



## LS-150™

Totally re-engineered for performance and accuracy, our premium-quality exam light delivers brilliant color and convenience.

- > Color-corrected 50-watt lamp delivers superior light quality for more accurate color rendition
- > Specially designed light head keeps luminaire cool to the touch, even during multiple exams
- > Small base profile frees up valuable exam room floor space
- > Encased, flexible 24 inch arm enables one-hand positioning—not possible with standard goose neck units
- > Multiple installation options to fit your needs

**44500** LS-150 Halogen Examination Light, Complete with Integral Floor Stand

**44510** LS-150 Halogen Examination Light, Wall/Table Mount Model

## LS-150™ Accessories

**44101** Table Clamp Accessory for 44510

**44350** Caster Base Accessory for 44500

**04450-U** 50 W Halogen Replacement Lamp

## Professional PenLite

For quick and effective evaluation of the pupillary reflex.

- > Bright, white halogen light for true tissue color rendition
- > Durable solid brass construction

**76600** Halogen Professional PenLite



**76600**  
Halogen Professional  
PenLite

## Professional PenLite Accessories

**07600-U** 2.5 V Halogen Lamp



**47300**  
EpiScope



**05232-U**  
Soft Carrying  
Case

## 3.5 V EpiScope® Skin Surface Microscope

Helps distinguish pigmented skin lesions such as basal cell carcinoma vs. malignant melanoma by providing the ability to magnify, illuminate and measure.

- > Provides 10 times magnification and illumination for crisp, clear viewing
- > Reticle contact plate facilitates measurement and tracking of lesion growth
- > Fits any 3.5 V Welch Allyn power source
- > Halogen illumination for true tissue color
- > Free reference guide available for patient education

**47300** EpiScope only

**47351** EpiScope Set with Rechargeable Handle

**47351-C** EpiScope Set with Convertible Handle

## 3.5 V EpiScope Skin Surface Microscope Accessories

**05232-U** Soft Carrying Case

**47320** Plano Contact Plate  
for EpiScope

**47310** Reticle Contact Plate  
for EpiScope

**03000-U** 3.5 V Halogen Lamp



## LS-200™

Bring high-quality, low-cost illumination to all your office-based procedures.

- > Superior maneuverability covers a wide range of positions without drifting or repositioning
- > 24 inch working distance with 2,500 footcandles light output and 3,200 °K color temperature enhances visualization and true tissue color
- > Three-lamp design ensures continuous, shadow-free illumination
- > 7 inch spot size with feathered edge definition and 12 inch depth of field reduce eye strain

**44200** LS-200 Light with Arm, Luminaire Section and Transformer (DOES NOT include mounting option)

## LS-200™ Accessories

**44250** LS-200 Mobile Stand

**44210** LS-200 Wall Mount

**06400-U** 25 W Halogen Lamp

# 140 Diagnostic Headlights

## LIGHTING



**46070**  
Solid-State Portable Headlight  
with Direct Power Supply/Charger

## Solid-State Portable Headlights

Provide illumination precisely where it's needed with cool, bright light that ensures comfort for both doctor and patient.

- > Weighing approximately four ounces, it's one of the most comfortable and lightweight headlights available
- > Fits inside a shirt pocket and allows complete freedom of movement
- > Supplies bright, white, shadow-free light for true tissue color
- > Solid-state lamp has 10,000-hour life expectancy with minimal degradation of light output

**46070** Solid-State Portable Headlight with Direct Power Supply/Charger

## Solid-State Portable Headlights Accessories

**05232-U** Soft Carrying Case

**49642** Soft Terrycloth Headband (case of 10)

**72270** Rechargeable Battery

**74170** Direct Power Source/Charger



Operate the Headlight from a Belt-Clip or Direct Power Source



**49020**  
Solid-State Procedure Headlight with Direct Power Supply

## Solid-State Procedure Headlights

Provide cool, bright white light that is ideal for a variety of office-based procedures. When you need targeted, well-defined illumination without the heat of traditional headlights, look to our solid-state headlights for superior quality and reliability.

- > Weighing approximately six ounces, it is extremely comfortable to wear for simple and complex procedures
- > Supplies bright, white, shadow-free light for true tissue color
- > Rechargeable battery pack makes this a completely portable, compact unit that allows freedom of movement
- > Enjoy the flexibility of operating the headlight using the direct power source or the rechargeable belt-clip portable power source
- > Solid-state lamp has 10,000-hour life expectancy with minimal degradation of light output

**49020** Solid-State Procedure Headlight with Direct Power Supply

## Solid-State Procedure Headlights Accessories

**49099** Carrying Case

**75260** Portable Power Source with Charger

**49000B** Procedure Headlight (no power supply)

**49097** Set, Front and Rear Pad

**72250** Replacement Cells for Battery Pack

**74180** Charger only for Portable Power Source

**75270** Direct Power Supply

**75201** Portable Power Source only



**49099**  
Carrying Case

## Solid-State Procedure Headlight Charger and Power Source



**75201**  
Portable Power Source only



**74180**  
Charger Only for Portable Power Source

## 142 Minor Procedure Headlights

### LIGHTING



**49520**  
Solarc MFI Headlight System  
(Light Source, Fiber Bundle,  
Headband/Luminaire)



**49500TW**  
Solarc MFI Light Source  
with Turret

### Solarc® MFI Headlight

Our durable, comfortable, and lightweight surgical headlight delivers the bright, white light expected from a xenon system, but at a fraction of the cost.

- > Welch Allyn's patented Solarc lamp provides white, xenon-quality light (5,500 °K) for more accurate color and clear tissue identification than 3,200 °K halogen systems
- > The Solarc arc lamp has a median life of 550 hours
- > Unique optical design provides at least 30% greater light intensity than leading halogen systems for excellent visualization of the operative site
- > Soft, flexible headband with a ratchet on the top and back for superior comfort and customized fit
- > Fiber optic bundle is encased in three layers of protective materials to eliminate fiber breakage and reduce costs

**49520** MFI Solarc® Headlight System (Light Source, Fiber Bundle, Headband/Luminaire)

**49520TW** MFI Solarc Headlight System with Turret (accepts Wolf, ACMI, Storz, Olympus)

### Solarc® MFI Headlight Accessories

**49500** Solarc MFI Light Source only

**49500TW** Solarc MFI Light Source with Turret (accepts Wolf, ACMI, Storz, Olympus)


**49550** Solarc MFI Mobile Stand

**49540** Solarc MFI Headlight Headband/Luminaire

**49541** Fiber Optic Bundle, 10 ft

**49543** Fiber Optic Bundle, 8 ft

**09500-U** Solarc MFI Replacement Lamp

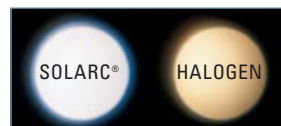

Welch Allyn Solarc lighting generates 50% whiter, brighter light for true tissue color