



DRYSTAR 5503

Imager

DRYSTAR 5503 IS A HIGH-THROUGHPUT, HIGH-RESOLUTION, MULTIPLE MEDIA SIZE DIRECT DIGITAL IMAGER THAT CAN BE EASILY NETWORKED FOR MAXIMUM PRODUCTIVITY.

- Outstanding imaging quality across all applications
- Multi-modality
- Multi-format with 3 media sizes on-line
- Diagnostic quality greyscale hardcopies of high quality
- Unique sorter function
- Mammography printing option



High resolution of 508 ppi

Improved workflow across departments

Offering high throughput of up to 100 14 x 17" sheets/hour, unique multi-format architecture and sorting by modality function, DRYSTAR 5503 is a versatile partner capable of handling a diversity of printing tasks from different sources with speed and reliability.

Every pixel tells a story

The Direct Digital Imaging technology employed by DRYSTAR 5503 is not only fast, clean and environmentally friendly, it is also one of the best ways of translating the high-resolution capability of 508 ppi at a spot size of 50 μ m to a final image. Agfa HealthCare's Direct Digital Imaging technology performs well beyond current industry standards, ensuring that every pixel in the image is fully software controlled for virtually flawless and consistent image quality. A significant reduction in wear and tear is due to an equally significant reduction in moving parts as a result of using Direct Digital Imaging technology.

Multi-format and 3 media sizes on-line

DRYSTAR 5503 features multiple-format handling, with the three most popular media sizes and/or types permanently on-line. The imager is thus capable of delivering CT, MRI, DSA, digital R&F, CR and DR images at high speed onto different DRYSTAR DT2 media. The DRYSTAR 5503 comes with three input trays. Every input tray can use five different film formats, ranging from 8 x 10" to 14 x 17", making the final image versatility of this stand-alone, small footprint unit remarkable. And what DRYSTAR 5503 gains in versatility, the user gains in convenience and time.



Three media trays for different on-line media sizes




Unique sorter eliminates high-traffic bottlenecks

Sorting it out

In order to increase efficiency and user-friendliness, DRYSTAR 5503 has a unique sorter function. This ensures that whichever modality is being served, whichever print command is currently being handled, all tasks will be carried out and sorted, according to modality or patient record, if available. With DRYSTAR 5503, high traffic density input no longer means bottlenecks and slow-speed output at the imager, virtually eliminating log jams.

Enhancing imaging quality

As part of our continual drive to provide you with perfect image quality, our award winning Direct Digital Imaging (DDI) technology has been enhanced to include A[#]Sharp technology. This technology intensifies imaging capability, enabling DDI to provide sharper image quality across all applications.

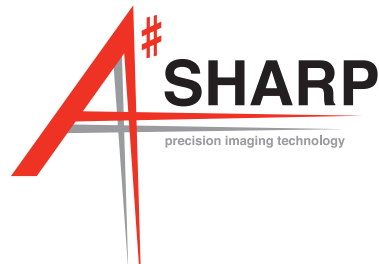

Integrated print solution

Through its intelligent matching of Direct Digital Imaging technology, media and imager, DRYSTAR 5503 is ideally suited to stand at the heart of an Agfa HealthCare integrated solution. Its state-of-the-art DRYSTAR DT2 media delivers superior, diagnostic quality greyscale hardcopies. Because it is heat sensitive rather than light sensitive, DT2 offers the convenience of daylight loading. The dry imaging technology means no more wet processing, no darkroom, no complicated adjustments or cleaning procedure, and no chemical disposal costs. DRYSTAR DT2 media can be used in all formats: 8 x 10", 10 x 12", 11 x 14", 14 x 14" and 14 x 17".



Increasing performance for mammo

As an option, DRYSTAR 5503 can also print mammography images. Using our dedicated DT2 Mammo film ensures you the highest image quality, in accordance with mammography standards. With its 508 dpi resolution and the A#Sharp technology, DRYSTAR 5503 produces the sharpest images, ideally suited for digital mammography. An automatic quality assurance tool is built in to help you with your daily quality checks. Images can be printed on 8 x 10" film, as well as on 10 x 12" film for modalities with larger detectors. As a third film size for mammography we also offer 11 x 14". In the mammography mode, the DRYSTAR 5503 can print up to 160 8 x 10" mammo films per hour, fitting easily into an environment with several modalities, while ensuring full performance for the mammography department.

By bringing together a diversity of modalities within an appropriate user-friendly integration – sharp high resolution imaging, state-of-the-art technology, excellent media and maximum user-friendliness – DRYSTAR 5503 is well equipped to serve the multiple demands of a busy department.

technical

SPECIFICATIONS

GENERAL

Dimensions & weight

- Dimensions: 72 x 71,5 x 141 cm (W x D x H)
(28,3 x 28,1 x 55,5 inch)
- Weight (without film): 193 kg (425 lb)

Power requirements

- Auto ranging 100 – 240 V: 50/60 Hz

Power consumption

- Average: 450 Watt
- Peak: 700 Watt
- Standby: 200 Watt

Media supply trays

- All three trays configurable for 100 sheets of
8 x 10 inch, 10 x 12 inch, 11 x 14 inch,
14 x 14 inch or 14 x 17 inch film

Operating conditions

- Temperature: 15 - 30 °C (59 - 86 °F)
- Humidity: 20 - 75% RH, non-condensing

Storage/Shipping conditions

- Temperature: -25 to +55 °C (-11 to 131 °F)
-40 to +70 °C for transport (-40 to 158 °F)
- Humidity: 10 - 95% RH, non-condensing

Heat dissipation

- Standby power: 200 Watt / 720 kJ/h
- Average printing power: 450 Watt / 1620 kJ/h
- Peak power: 700 Watt / 2520 kJ/h

Safety

- IEC 60601-1 + A1 + A2
- EN 60601-1 + A1 + A2
- UL 60601-1
- CSA 22.2 no. 601.1-M90
- GB 4943-2001

PERFORMANCE

Throughput

- 8 x 10 inch: 160 sheets/hour (23 sec. per sheet)
- 14 x 17 inch: 100 sheets/hour (36 sec. per sheet)
- Access time first sheet: 57 sec. (8 x 10 inch)
76 sec. (14 x 17 inch)

Addressable print area

- 8 x 10 inch: 3892 x 4920 pixels
- 14 x 17 inch: 6962 x 8408 pixels

Printing resolution

- Geometrical: 508 ppi
- Spot size: 50 µm
- Contrast: 14 bit printing

Connectivity

- Ethernet TCP/IP
- Protocol: DICOM 3.0

MEDIA

Types

- DRYSTAR DT2 B: blue base
- DRYSTAR DT2 C: clear base
- DRYSTAR DT2 Mammo

Sizes

- 8 x 10 inch, 10 x 12 inch, 11 x 14 inch, 14 x 14 inch
and 14 x 17 inch for DRYSTAR DT2 B/C
- 8 x 10 inch, 10 x 12 inch,
11 x 14 inch for DT2 Mammo
- Three on-line sizes

Agfa and the Agfa rhombus are trademarks of Agfa-Gevaert N.V., Belgium, or its affiliates. DRYSTAR is a trademark of Agfa HealthCare NV, Belgium or its affiliates. All other trademarks are held by their respective owners and are used in an editorial fashion with no intention of infringement. The data in this publication are for illustration purposes only and do not necessarily represent standards or specifications, which must be met by Agfa HealthCare. All information contained herein is intended for guidance purposes only, and characteristics of the products and services described in this publication can be changed at any time without notice. Products and services may not be available for your local area. Please contact your local sales representative at agfa.com for availability information. Agfa HealthCare diligently strives to provide as accurate information as possible, but shall not be responsible for any typographical error.

© 2010 Agfa HealthCare NV

All rights reserved

Printed in Belgium

Published by Agfa HealthCare NV

B-2640 Mortsel - Belgium

5PQNA US 00201004