

# Sterile Fields, Towels and Drapes



## Sterile Fields

18" x 26" sterile. Tissue/poly/tissue. Blue/white.


**CATALOG NO. 696** 

*Packed 300/carton (6 dispensers of 50 ea.)*

18" x 26" non-sterile. Tissue/poly/tissue. Autoclavable: polypropylene film center. For use with individual treatment trays. Blue/white.

**CATALOG NO. 695** 

*Bulk packed, 500/carton*

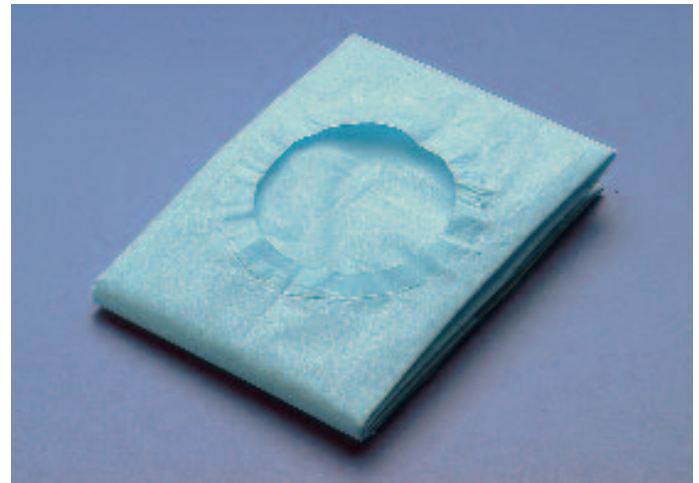

## Minor Procedure Fenestrated Drapes

18" x 26" sterile. Tissue/poly/tissue, with 3" round fenestration. Blue/white.

**CATALOG NO. 697** 

*Packed 300/carton (6 dispensers of 50 ea.)*

18" x 26" non-sterile. Tissue/poly/tissue with 3" round fenestration. Autoclavable: polypropylene film center. For use with individual treatment trays. White.

**CATALOG NO. 694** 

*Bulk packed, 300/carton*

We offer a variety of towels and drapes to meet the many needs of hospitals and clinics.

Our versatile "tissue/poly/tissue" towels and drapes feature a polyethylene lining between two layers of absorbent strength tissue, virtually eliminating the risk of bacterial migration and strike through. They create a convenient sterile field, and, for minor surgical procedures, are available with a 3" fenestration, with or without an adhesive patch around the fenestration.


## HERE ARE SOME OF THE MANY FEATURES YOU'LL FIND IN OUR STERILE FIELDS:



Polylining, affixed by a heat sealing process, prevents passage of liquids, air, bacteria – preserves integrity of sterile field.



All sterile fields are available in convenient dispenser boxes of 50 towels each.

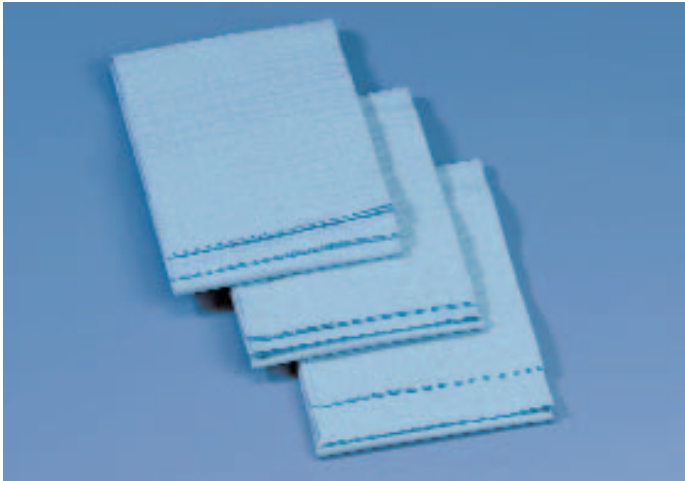


### Kaycel® Towels – Extra Strength

19½" x 28" non-sterile, absorbent, blue general purpose Kaycel® towel.

**CATALOG NO. 6935** 

*Bulk packed 500/carton*

19½" x 23" non-sterile, absorbent, general purpose blue Kaycel® towel.

**CATALOG NO. 6710** 

*Bulk packed 200/carton*

Absorbent and non-absorbent Kaycel® drapes and towels are available with or without an adhesive strip.

Our non-sterile towels and drapes with polypropylene linings are completely autoclavable, making them ideal for hospital-prepared trays.




...WHEN AN **EMBOSS**ED, HIGHER LOFT DRAPE IS DESIRED

### Non-Sterile Embossed Drape

**CATALOG NO. 3270**  Blue

*Bulk packed, 300 per carton*

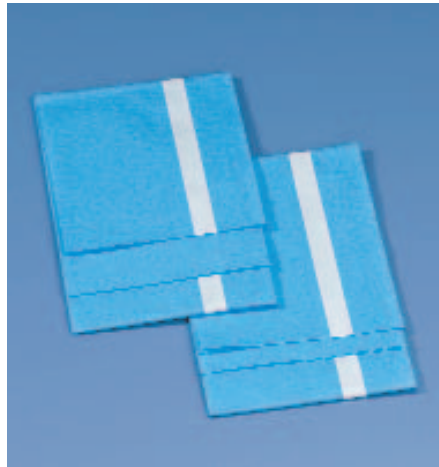

### Utility Towels and Drapes

(2)18" x 26" sterile, non-absorbent, blue utility drapes with adhesive strip.


**CATALOG NO. 680** 

*Packed 2 drapes/package (50 packages/carton)*

18" x 26" sterile. Tissue/poly/tissue, with adhesive patch around 3" fenestration. White.

**CATALOG NO. 698** 

*Packed 200/carton (4 dispensers of 50 ea.)*

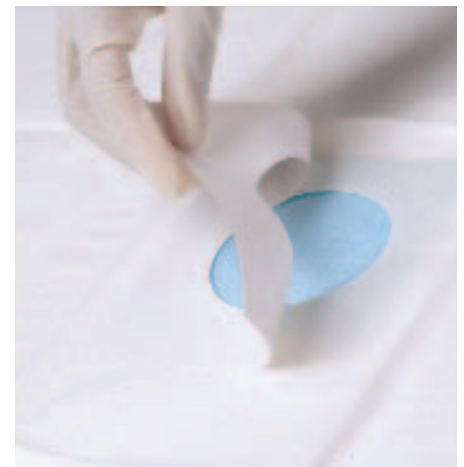

Fenestration is die-cut and heat sealed – leaves a smooth, even edge – prevents linting. Shown above, adhesive patch around fenestration of our tissue/poly/tissue drape.