



**Think it's easy to comply with infection control guidelines at your facility? Think again.**

Even the most scrupulously maintained healthcare facility can be teeming with germs that can cause healthcare-associated infections in patients, staff, and visitors.

Outside contaminants that are brought into the environment on a daily basis and contaminants from the patients themselves are a constant hazard. Something as seemingly innocent as a visitor's hug, or even a clinician's healing touch, can leave behind dangerous germs.


To help prevent the spread of healthcare-associated infections, the CDC recommends that clinicians sanitize their hands and don fresh personal protective equipment (PPE) before seeing each patient, in the following sequence<sup>1</sup>:

- Step 1:** Sanitize hands
- Step 2:** Don gown
- Step 3:** Don mask
- Step 4:** Don gloves

**THE KIMBERLY-CLARK\* PPE Dispensing System makes it easier for your staff to comply!**



*Under special lighting, with a luminescent substance to represent bacteria, you can see how easy it is to transfer germs from clinician to patient with a single touch.*

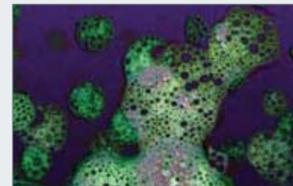

*This sequence shows how a facemask, after a single use, can become thoroughly contaminated with germs.*

*At Kimberly-Clark, our mission is to deliver clinical solutions that you can depend on to meet the demands of your fast-paced world. Whether your needs involve preventing healthcare-associated infections, surgical and digestive solutions or pain management, with Kimberly-Clark you'll always have one less worry.*



Healthcare-Associated Infection Solutions



Surgical Solutions



Digestive Health



Pain Management



**Commitment to Excellence**

If, for any reason, our products do not meet your expectations, please let us know your comments or suggestions for improvement. Your input will result in a concerted effort on our part to meet your requirements. Our goal is to provide quality products that completely meet your needs time after time.

For more information, please call your Kimberly-Clark representative or visit our web site at [www.kchealthcare.com](http://www.kchealthcare.com).

\* Registered Trademark or Trademark of Kimberly-Clark Worldwide, Inc. The COLOR PURPLE is a Registered Trademark of Kimberly-Clark Worldwide, Inc. ©2007 KCWW. All rights reserved. H00187

1. Guidance for the Selection and Use of Personal Protective Equipment (PPE) in Healthcare Settings. CDC, Dec 2, 2005



**The KIMBERLY-CLARK\* PPE Dispensing System:**  
*One less worry about infection control*

# The KIMBERLY-CLARK® PPE Dispensing System makes it easy to comply with the CDC Infection Control Guidelines.

**KIMBERLY-CLARK®**  
PPE Dispensing System

The KIMBERLY-CLARK® PPE Dispensing System is an all-in-one, total solution for dispensing and storing personal protective equipment and hand sanitizer; products you know and trust, only from Kimberly-Clark.

- Provides convenient access to PPE
- Promotes compliant infection control practices
- Helps reduce the spread of HAIs

## It's modular...it's portable...and it makes Kimberly-Clark PPE accessible throughout your facility.

The KIMBERLY-CLARK® PPE Dispensing System promotes compliant infection control practices to help reduce HAIs. The system makes it simple for clinical staff to sanitize their hands and change PPE between patients.

The dispensing module can be personalized with your facility's name and can be wall-mounted or attached to a wheeled storage module, to travel with staff wherever it's needed.



## Dispensing

### The all-in-one, total PPE solution...

The KIMBERLY-CLARK® PPE Dispensing System's ergonomically designed dispensing module features easy-to-follow icons that promote hand sanitization and proper donning techniques. And it's designed to accommodate the Kimberly-Clark PPE that staff members use most:

- KMCARE ANTIBACTERIAL™ ALL-N-1 Hand Sanitizer
- CONTROL® Cover Gowns
- STERLING® Nitrile Exam Gloves and/or PURPLE NITRILE® Gloves
- KIMBERLY-CLARK® Procedure Masks
- FLUIDSHIELD® Visor Masks

### ...that helps reduce contamination.

It dispenses PPE using a unique, proprietary system that helps prevent contamination during the selection and donning process.

**Gloves and procedure masks are dispensed using a spring-assisted mechanism.** When combined with Kimberly-Clark's proprietary box openings, the need to reach into the boxes is eliminated. This helps minimize the risk of contamination that can occur when normal dispensing requires a user's hand to go deep into a box to retrieve gloves or masks.

**FLUIDSHIELD® Visor Masks are dispensed mask first,** so that visors stay virtually fingerprint-free.

**Exam gloves are dispensed cuff first,** which makes donning faster and easier...as well as helping to prevent contamination of the glove surface.



The unit can be personalized with your facility's name.

Easy to follow icons promote proper donning techniques.

### Dispensing Module

The dispensing module can be wall mounted or attached to the wheeled storage module.

### Storage Module

## Storage

The wheeled storage module can be used to hold extra quantities of Kimberly-Clark items for rapid refill. Additional boxes of gloves and masks are also conveniently stored in the dispensing module, behind the box currently being dispensed.

### Loading the dispensing module is simple.

The glove and mask boxes featuring spring-assisted dispensing open on both sides. You simply position the open back of the box in front of the spring and close the unit. Gowns lie flat in the dispensing module. And the hand sanitizer container fits easily onto the bracket mount on the dispensing unit.



Exam gloves dispense cuff first, preventing contamination of the glove surface.



Spring-assisted dispensing minimizes the risk of contamination.



The wheeled storage module provides plenty of extra room for additional supplies.

For more information about the KIMBERLY-CLARK® PPE Dispensing System, contact your Kimberly-Clark Sales Representative, call 1-800-KCHELPS, or visit our website at [www.kchealthcare.com](http://www.kchealthcare.com).

The KIMBERLY-CLARK® PPE Dispensing System is just one of the clinical solutions you can depend on to meet the demands of your fast-paced world.

Whether your needs involve preventing healthcare-associated infections, surgical and digestive solutions, or pain management, with Kimberly-Clark, you'll always have one less worry.

## Product Details

Dimensions:	w	x	d	x	h
Storage Module:	30"	x	20"	x	26.3"
Dispensing Module:	19.6"	x	17.3"	x	50.5"
Total System:	30"	x	20"	x	76.1"

## Capacity:

### Dispensing Module:

- CONTROL® Gowns - 30 each (universal)
- STERLING® Nitrile/PURPLE NITRILE® Gloves - 6 BOXES (3 dispensing, 3 in reserve)
- FLUIDSHIELD® Visor Masks - 1 box<sup>†</sup>  
Procedure Masks - 1 box<sup>†</sup>
- KMCARE ANTIBACTERIAL™ All-N-1 Hand Sanitizer - 1- 1.2L Container

<sup>†</sup> One box of either type of mask can be stored in reserve.

### Storage Module:

- Holds additional products for restocking Dispensing Module
- Flexible, to accommodate products you use the most

## Product Codes:

System	Code	Description
29901	Customized Dispensing Module	
29902	Branded Dispensing Module	
29903	Storage Module	

## Gowns

69979	CONTROL® Yellow, Universal
69981	CONTROL® Blue, Universal
69983	CONTROL® White, Universal

## Gloves

STERLING® Nitrile for use with the PPE Dispensing System	
55095	Small
55096	Medium
55097	Large

PURPLE NITRILE® for use with the PPE Dispensing System	
55101	Small
55102	Medium
55103	Large

## Masks

47149	FLUIDSHIELD® with Visor for use with the PPE Dispensing System
47119	Procedure Mask for use with the PPE Dispensing System

## Hand Sanitizer

95150	Instant Hand Sanitizer
95151	Moisturizing Instant Hand Sanitizer