



Proxima

Surgical Gowns

Designed, Tested and Trusted

by O.R. Professionals

Proxima Surgical Gowns

Designed, Tested and Trusted by O.R. Professionals

Proxima was designed by you.

A surgical gown is the uniform of the surgical team. Just as athletes depend on their apparel to enhance their performance, you should expect the same from your surgical gown.

It should offer you uncompromising protection. It should be comfortable. And it should be the last thing you have to worry about during a surgical procedure.

When we developed our Proxima line of surgical gowns we went straight to the people who use them. We spoke to O.R. professionals across the country to get a sense of what features they liked and disliked about other surgical gowns, and then worked to develop the best gowns.

The power of choice.

We recognize that your O.R. team may have different preferences in a surgical gown. That's why we offer several styles of gowns. Each brand of gown features a different fabric with its own unique qualities. However, there are a set of benefits that all these gowns share.



Proxima gowns feature a velcro-style neck closure for easy adjustability.

All Proxima gowns are cut wider in the chest and sleeves for a better fit. Instead of using snaps, our gowns feature an adjustable velcro-style neck closure. This allows the neckline to easily adjust to a wide vari-

ety of people, while also helping to prevent gapping in the front which serves as additional protection against fluids.

All Proxima gowns meet the Standard for Flammability of Clothing Textiles CPSC 16 CFR Part 1610. In addition, all gowns meet and/or exceed the latest recommendations from AAMI (Association for the Advancement of Medical Instrumentation).

No other manufacturer offers this kind of choice.

All Proxima gowns are latex free.

Keeping you and the patient in mind, every gown is latex free. All of our packaging is clearly labeled with our latex-free symbol for added peace of mind.



For more information, contact your Medline representative or call us at 1-800-MEDLINE.



Proxima

Aurora

Made with Suprel® medical fabric from DuPont™

Suprel medical fabric is based on DuPont's new Nonwovens Advanced Composite technology. This bi-component fabric is made of polyester (for strength) and polyethylene (for comfort).



The innovative polyester and polyethylene combination as viewed under an electron microscope.

Uncompromising performance

Aurora gowns offer a high level of protection without compromising comfort. Featuring new fabric technology, Aurora raises the bar for single-use O.R. gowns.

Breathable and comfortable

Aurora gowns have a smooth texture that glides easily without catching or grabbing, giving you incredible freedom of movement. Plus, the gown helps to quickly transfer heat away from your body, keeping you cool and comfortable.

Excellent protection

For a high level of protection, Aurora gowns are a great choice. The gown's high hydrostatic head offers superior barrier performance.



Aurora™ Surgical Gowns


Aurora gowns offer an unbeatable combination of protection and comfort.



Ordering Information

Non-Reinforced

DYNJP2701	Large (43", 109 cm), Non-Sterile Item #: SPT-2701	30/cs
DYNJP2702	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2702	30/cs
DYNJP2703	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2703	18/cs
DYNJP2715	Small (41", 104 cm), Non-Sterile Item #: SPT-2715	30/cs

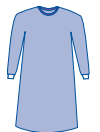

AAMI Level 3

Non-Reinforced, Raglan Sleeve


DYNJP2801	Large (45", 114 cm), Non-Sterile Item #: SPT-2801	32/cs
DYNJP2802	X-Large (48", 122 cm), Non-Sterile Item #: SPT-2802	30/cs
DYNJP2803	XX-Large (50", 127 cm), Non-Sterile Item #: SPT-2803	28/cs



AAMI Level 3

Fabric-Reinforced

DYNJP2704	Large (43", 109 cm), Non-Sterile Item #: SPT-2704	30/cs
DYNJP2705	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2705	30/cs
DYNJP2706	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2706	18/cs

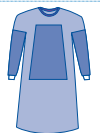

AAMI Level 3

Fabric-Reinforced, Raglan Sleeve

DYNJP2804	Large (45", 114 cm), Non-Sterile Item #: SPT-2804	30/cs
DYNJP2805	X-Large (48", 122 cm), Non-Sterile Item #: SPT-2805	28/cs
DYNJP2806	XX-Large (50", 127 cm), Non-Sterile Item #: SPT-2806	28/cs



AAMI Level 3

Poly-Reinforced

DYNJP2707	Large (43", 109 cm), Non-Sterile Item #: SPT-2707	30/cs
DYNJP2708	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2708	30/cs
DYNJP2709	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2709	18/cs



AAMI Level 4

Poly-Reinforced, Extra-Long

DYNJP2724	Large (52", 132 cm), Non-Sterile Item #: SPT-2724	28/cs
DYNJP2725	X-Large (56", 142 cm), Non-Sterile Item #: SPT-2725	28/cs
DYNJP2726	XX-Large (56", 142 cm), Non-Sterile Item #: SPT-2726	18/cs



AAMI Level 4

Poly-Reinforced, Raglan Sleeve

DYNJP2807	Large (45", 114 cm), Non-Sterile Item #: SPT-2807	32/cs
DYNJP2808	X-Large (48", 122 cm), Non-Sterile Item #: SPT-2808	30/cs
DYNJP2809	XX-Large (50", 127 cm), Non-Sterile Item #: SPT-2809	28/cs



AAMI Level 4

Poly-Reinforced, Breathable Impervious Sleeves

DYNJP2727	Large (43", 109 cm), Non-Sterile Item #: SPT-2727	30/cs
DYNJP2728	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2728	30/cs
DYNJP2729	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2729	18/cs
DYNJP2712	XXX-Large (49", 124 cm), Non-Sterile Item #: SPT-2712	18/cs



AAMI Level 4



Eclipse

Made with breathable Softesse® medical fabric from DuPont™

For more than 25 years, DuPont Softesse medical fabric has been considered an industry standard for medical gowns and drapes. The fabric features advanced material technology that provides a splash-resistant barrier that's breathable. It keeps you cool and comfortable while you work, helping you feel confident and in control.

Comfort you can feel

The first feature you will notice about the Eclipse gown is the way it feels. The soft, cloth-like texture says "comfort" as soon as you put it on. Those who spend much of their day in a surgical gown will appreciate the comfort of Eclipse.

Variety of styles

Eclipse gowns feature a variety of styles and protection levels including non-reinforced, fabric-reinforced, Cool-Zone and poly reinforced gowns that are also available with breathable sleeves and extra-long sizing.



Eclipse[®] Surgical Gowns

Eclipse gowns offer everyday protection and comfort.



Ordering Information

Non-Reinforced

DYNJP2001	Large (43", 109 cm), Non-Sterile Item #: SPT-2001	30/cs
DYNJP2002	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2002	30/cs
DYNJP2003	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2003	18/cs
DYNJP2004	XXX-Large (50", 127 cm), Non-Sterile Item #: SPT-2004	18/cs
DYNJP2005	Small (41", 104 cm), Non-Sterile Item #: SPT-2005	30/cs
DYNJP2009	XXXX-Large (49", 124 cm), Not available non-sterile	18/cs



AAMI Level 2

Fabric-Reinforced

DYNJP2101	Large (43", 109 cm), Non-Sterile Item #: SPT-2101	30/cs
DYNJP2102	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2102	30/cs
DYNJP2103	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2103	18/cs



AAMI Level 2

Cool-Zone, Fabric-Reinforced

DYNJP2104	Large (43", 109 cm), Non-Sterile Item #: SPT-2104	30/cs
DYNJP2105	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2105	30/cs

Poly-Reinforced

DYNJP2201	Large (43", 109 cm), Non-Sterile Item #: SPT-2201	30/cs
DYNJP2202	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2202	30/cs
DYNJP2203	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2203	18/cs

Cool-Zone, Poly-Reinforced

DYNJP2213	Large (43", 109 cm), Non-Sterile Item #: SPT-2213	30/cs
DYNJP2214	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2214	30/cs



AAMI Level 4

Poly-Reinforced, Extra-Long

DYNJP2204	Large, Extra-Long (52", 132 cm), Non-Sterile Item #: SPT-2204	28/cs
DYNJP2205	X-Large, Extra-Long (56", 142 cm), Non-Sterile Item #: SPT-2205	28/cs
DYNJP2206	XX-Large, Extra-Long (56", 142 cm), Non-Sterile Item #: SPT-2206	18/cs

Poly-Reinforced, Breathable Impervious Sleeve

DYNJP2207	Large (43", 109 cm), Non-Sterile Item #: SPT-2207	30/cs
DYNJP2208	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2208	30/cs
DYNJP2209	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2209	18/cs
DYNJP2212	XXX-Large (50", 127 cm), Non-Sterile Item #: SPT-2212	18/cs



AAMI Level 4



Prevention Plus

Advanced technology provides best combination of protection and comfort

Prevention Plus gowns offer a high level of protection while providing increased breathability and comfort. The gown is comprised of three layers: a patent-pending hybrid film laminated between two layers of bi-component spunbond material. Because moisture vapor is thousands of times smaller than the smallest known virus, the microporous film allows moisture vapor to escape while preventing transmission of fluids and increasing comfort.

Impervious protection

Prevention Plus material meets ASTM test for viral penetration (ASTM F1671) and ASTM test for blood penetration (ASTM F1670).

Generous Sizing

Prevention Plus gowns are cut wider through the chest and sleeves for a roomier fit. Plus, they are available in both extra-long and standard lengths to accommodate a variety of needs in the O.R.



Prevention™ Plus Surgical Gowns

Prevention Plus gowns offer the highest level of protection with breathable microporous technology.

Ordering Information

DYNJP2301P	Large (45", 114 cm), Non-Sterile Item #: 19310P	24/cs
DYNJP2302P	X-Large, Extra-Long (56", 142 cm) Non-Sterile Item #: 21394P	24/cs
DYNJP2303P	XX-Large, Extra-Long (56", 142 cm), Non-Sterile Item #: SPT-2303P	22/cs
DYNJP2306P	Large (43", 109 cm) Non-Sterile Item #: SPT-2306P	24/cs
DYNJP2307P	X-Large (47", 119 cm) Non-Sterile Item #: SPT-2307P	24/cs
DYNJP2308P	XX-Large (49", 124 cm) Non-Sterile Item #: SPT-2308P	22/cs




AAMI Level 4

Proxima

Sirus

Four layers of protection

Our SMS technology utilizes four layers of fabric to provide a high level of fluid repellency. This helps provide a barrier to blood, body fluids and pathogens. The fabric features two densely-packed meltblown layers sandwiched between two strong, spunbond outer layers.

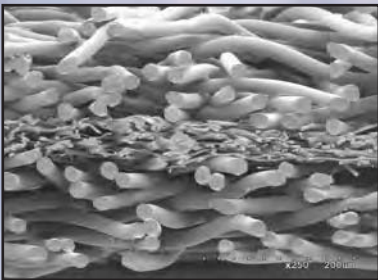

SMS layers as viewed under a microscope



Sirus™ Surgical Gowns

Sirus gowns offer a level of comfort superior to other SMS gowns.

Ordering Information

Non-Reinforced

DYNJP2001S	Large (43", 109 cm), Non-Sterile Item #: SPT-2001SMS	20/cs
DYNJP2002S	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2002SMS	20/cs
DYNJP2003S	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2003SMS	18/cs



AAMI Level 3

Non-Reinforced, Raglan Sleeve

DYNJP2401	Large (45", 114 cm), Non-Sterile Item #: SPT-2401SMS	32/cs
DYNJP2402	X-Large (48", 122 cm), Non-Sterile Item #: SPT-2402SMS	30/cs
DYNJP2403	XX-Large (50", 127 cm), Non-Sterile Item #: SPT-2403SMS	28/cs



AAMI Level 3

Fabric-Reinforced

DYNJP2101S	Large (43", 109 cm), Non-Sterile Item #: SPT-2101SMS	20/cs
DYNJP2102S	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2102SMS	20/cs
DYNJP2103S	XX-Large (49", 124 cm), Non-Sterile Item #: SPT-2103SMS	18/cs



AAMI Level 3

Fabric-Reinforced, Raglan Sleeve

DYNJP2501	Large (45", 114 cm), Non-Sterile Item #: SPT-2501SMS	30/cs
DYNJP2502	X-Large (48", 122 cm), Non-Sterile Item #: SPT-2502SMS	28/cs
DYNJP2503	XX-Large (50", 127 cm), Non-Sterile Item #: SPT-2503SMS	28/cs



AAMI Level 3

Poly-Reinforced

DYNJP2201S	Large (43", 109 cm), Non-Sterile Item #: SPT-2101SMS	20/cs
DYNJP2202S	X-Large (47", 119 cm), Non-Sterile Item #: SPT-2102SMS	20/cs
DYNJP2203S	X-Large (49", 124 cm), Non-Sterile Item #: SPT-2103SMS	18/cs



AAMI Level 4

Poly-Reinforced, Extra-Long

DYNJP2204S	Large, Extra-Long (52", 132 cm), Non-Sterile Item #: SPT-2204SMS	20/cs
DYNJP2205S	X-Large, Extra-Long (56", 142 cm), Non-Sterile Item #: SPT-2205SMS	20/cs
DYNJP2206S	XX-Large, Extra-Long (56", 142 cm), Non-Sterile Item #: SPT-2206SMS	18/cs

Poly-Reinforced, Raglan Sleeve

DYNJP2601	Large (45", 114 cm), Non-Sterile Item #: SPT-2601SMS	32/cs
DYNJP2602	X-Large (48", 122 cm), Non-Sterile Item #: SPT-2602SMS	30/cs
DYNJP2603	XX-Large (50", 127 cm), Non-Sterile Item #: SPT-2603SMS	28/cs



AAMI Level 4