

Mobile Adjustable Height Clear Lead Acrylic Barrier



The Adjustable Height Barrier allows the user to position the clear pb at any height level between 45 1/2" and 68". By use of an internal counter-balancing mechanism, the clear pb can be moved up or down with minimal effort.

The overall width of the base is 31 inches and provides 1mm pb equivalent protection. The upper 28 inch wide clear pb provides 0.5mm lead equivalency protection.

This strong yet lightweight barrier is convenient to use for dozens of applications. The 7" legs and casters facilitate easy transport ensuring you will get maximum usage from it. When completely closed the barrier is easy to store. **#80-610.**

Mobile Non-Adjustable Clear Acrylic X-ray Barriers



Wolf offers 4 styles of our Non-adjustable Mobile Barriers, (see the table below for the sizes and lead equivalency protections available). Designed with a shatter resistant clear pb viewing area at the top and a white opaque bottom that is also lead lined for your complete protection, the framework is stainless steel with a brushed finish. Casters lock in place for stable positioning.

	Regular 80-602	Jumbo 80-606	Deluxe 80-600	Deluxe 1.0mm 80-605
Overall Size	75"h x 30"w	75"h x 30"w	75"h x 48"w	75"h x 30"w
Window Panel Lead Equiv.	24"h x 30"w 0.5mm	48"h x 30"w 0.5mm	30"h x 48"w 0.5mm	36"h x 48"w 0.5mm
Opaque Panel Lead Equiv.	48"h x 30"w 0.8mm	24"h x 30"w 0.8mm	42"h x 48"w 0.8mm	36"h x 48"w 0.8mm

Adjustable Mobile Shield




Designed to protect the patient against scattered radiation, the Wolf Mobile Shield is ideal for chest, skull and spine radiography. The screen is made of lead vinyl with a .5mm lead equivalency protection and is attractively covered on both sides by easy to maintain vinyl. It is 24" (60cm) wide and is adjustable with a minimum height of 24" (60cm) and a maximum height of 48" (122cm). It's mounted on a durable white finished frame and rolls easily to any location on 4 conductive casters. Specifications: 48" H (open) x 30" H (closed) x 26 3/4" W x 18"D (at base). **#23111**

Clear Pb Mobile Radiation Shield



The Clear Pb Mobile X-Ray Barrier provides a full-body shield for the doctor or R.T. The optically clear, lead-plastic viewing panel provides exceptional visibility and 0.5 mm lead-equivalent protection. The shatter-resistance and durability of Clear-Pb ensures longlife under heavy use. The Clear Pb shield can safely shield one or two persons. The Clear Pb window rests on an opaque lower panel, casters allow the barrier to be moved easily. Specifications: 75"H x 30"W x 24" D. Weight: 113 lbs.

#56604. (Please specify the lead equivalency desired when ordering.)

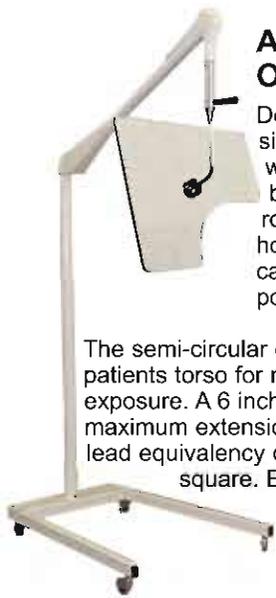

Adjustable Over Table Barrier

Designed for protection in any situation, our Over-Table Barrier wheels on any location. The barrier mount allows 360° rotation on vertical and horizontal planes so the barrier can be used in unlimited working positions.

The semi-circular cutout fits closely to the patients torso for maximum shielding from x-ray exposure. A 6 inch offset pivot provides a maximum extension of 36". The barrier is 0.5mm lead equivalency clear acrylic and is 24 inches square. Barrier height is adjustable from 63" above floor to 30" at the center. **#80-615**



Center Mount Barrier

Our "CMB's" have a ball and socket center mount that rotates 360°. It is supplied with your choice of a single or dual mount ceiling column or can be used with our ceiling track shown above right. The barrier is 76cm high and 61cm wide and provides 0.5mm lead equivalency protection. **#80-596.**



WolfLite Lamp With Arm Mounting

The WolfLite lamp runs on AC, 110-120V power and features adjustable focus with a light intensity of 35K Lux from one 24 Volt, 50 watt Halogen bulb. The lamp head is 22cm in diameter. Focusable illuminated field is 14 to 25 cm. Can be used in combination with any ceiling mount barrier. **#80-511.**

To order a custom shield of any type fill out the form on pg. 50

Ceiling Track System



Our Ceiling Track is made to be used with our "CMB" Mounted Barrier (80-596) or our Overhead Barrier (80-568). Shown below.

The perfect accessory for getting maximum usage from your barrier system, the track allows you to slide your barrier more than 8 feet. The track is 2.5 meters long (98.5") and comes equipped with one carriage that handles one column. **#80-270**



Frame Mounted Overhead Barriers

Can be used with a single or a dual ceiling mount, or attached to the ceiling track shown above. All options allow 360° barrier rotation. 0.5mm lead equivalent.

Description	Dimensions (HxW)	Order #
With cutout	23.5 x 31.5"	80-568
Dual Ceiling Mount	15.75" x 19.5"	80-645



Protective Curtains For Overhead Barriers

These curtains (shown attached to the bottom of a barrier) add an extra 8" of protection at the bottom of the barrier, easy to install via supplied velcro.

For Barrier **#80-654**.....Curtain **#80-356**
 For Barrier **#80-596**.....Curtain **#80-361**
 For Barrier **#80-568**.....Curtain **#80-380C**



3 Section Double Pivot Under Table Shield

This 32"W x 28" L shield comes with a removable 9" high top. It's equipped with a bracket that locks onto a table side accessory rail up to 28.5 H x 8.8mm D. The flexible 0.5mm lead equivalent overlapping curtains are attached to a pivot arm. **#80-404.**



Single Section Single Pivot Shield

24"W x 28"H with a removable 9" high top shield the full 24" length. It has a single pivoting arm with 0.5mm lead equ. protection. Equipped with a 17cm wide mounting bracket that locks into a table side accessory rail up to 28.5mm high by 8.5mm deep. **#80-403**